


Introducing the new **Dual** CHEMRxnHUB

YOUR **Rx** FOR PROCESS CHEMISTRY

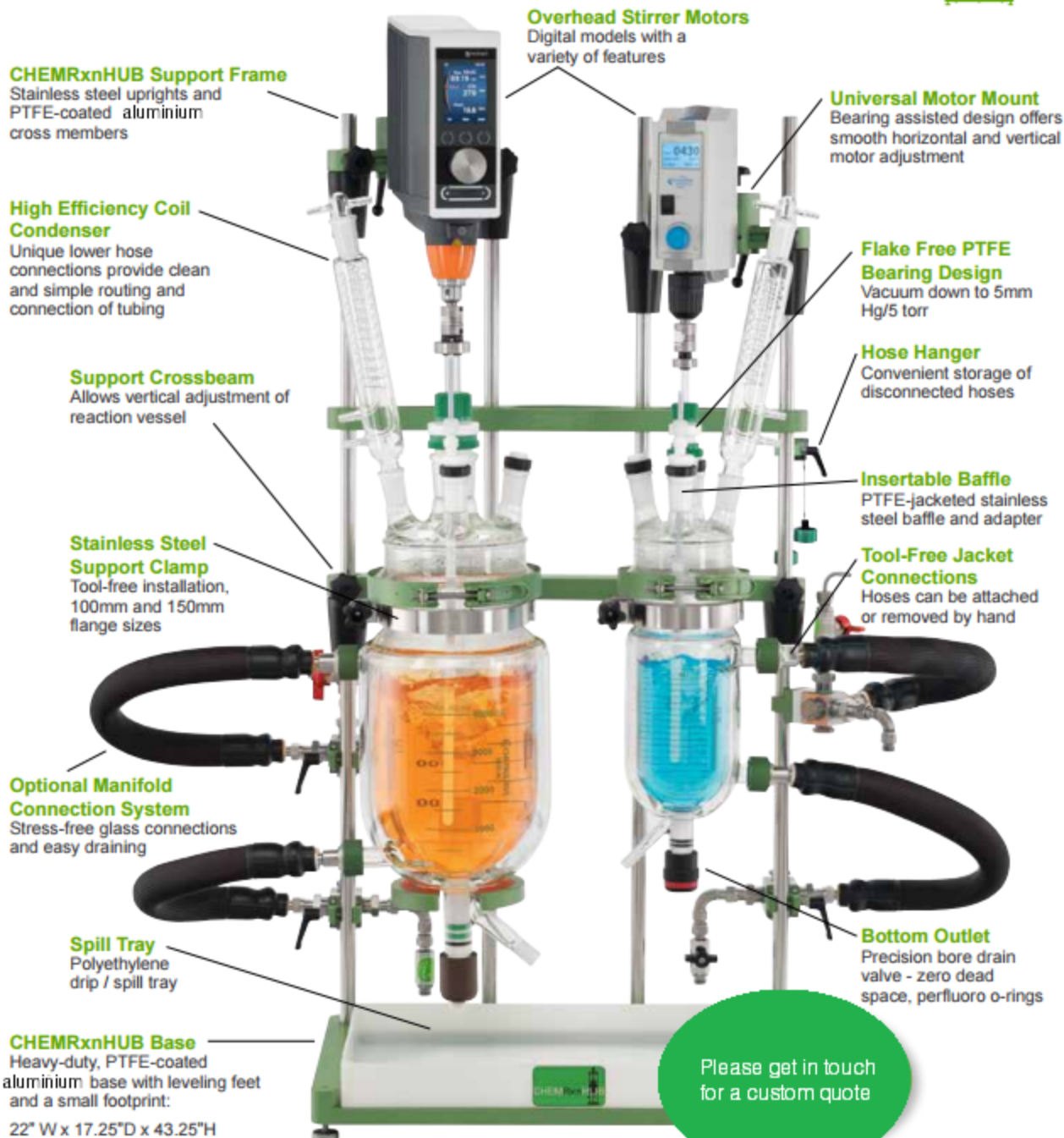

CHEMRxnHUB Support Frame

Stainless steel uprights and PTFE-coated aluminium cross members

High Efficiency Coil Condenser

Unique lower hose connections provide clean and simple routing and connection of tubing

Support Crossbeam

Allows vertical adjustment of reaction vessel

Stainless Steel Support Clamp

Tool-free installation, 100mm and 150mm flange sizes

Optional Manifold Connection System

Stress-free glass connections and easy draining

Spill Tray

Polyethylene drip / spill tray

CHEMRxnHUB Base

Heavy-duty, PTFE-coated aluminium base with leveling feet and a small footprint:

22" W x 17.25"D x 43.25"H

Overhead Stirrer Motors

Digital models with a variety of features

Universal Motor Mount

Bearing assisted design offers smooth horizontal and vertical motor adjustment

Flake Free PTFE Bearing Design

Vacuum down to 5mm Hg/5 torr

Hose Hanger

Convenient storage of disconnected hoses

Insertable Baffle

PTFE-jacketed stainless steel baffle and adapter

Tool-Free Jacket Connections

Hoses can be attached or removed by hand

Bottom Outlet

Precision bore drain valve - zero dead space, perfluoro o-rings

Please get in touch for a custom quote