


# enbio

To whom it may concern,

Due to the innovative and unique DET technology of the Microjet microwave autoclave, we recommend adhering to the sterilisation volumes mentioned in the table below.

**TABLE 1. RECOMMENDED STERILISATION VOLUMES**

Container - nominal capacity (ml)	Minimal volume (ml)	Maximum volume (ml)
100	100	100
200	100	150
250	150	200
500	150	400
750	250	500

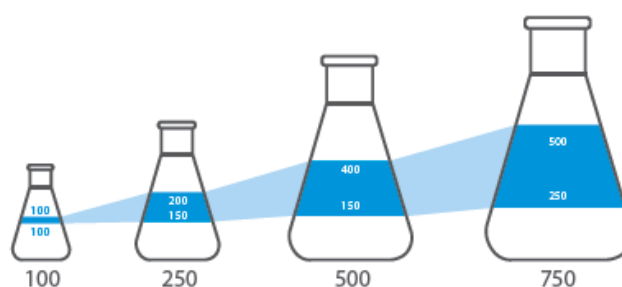

We recommend using borosilicate glass i.e DURAN, SIMAX bottles with a blue screw [heat resistant up to 140°C] and narrow neck Erlenmeyer flasks which are covered with aluminium foil.

Volume Range (ml)	Bottle capacity (ml)	Erlenmeyer flask capacity (ml)	DURAN cat no (bottle/flask)
100	100	100	21801245 / 2121624
100-150	200	200	flask only 2121632
150-200	250	250	21801365 / 2121636
150-400	500	500	21801445 / 2121644
250-500	750	750	bottle only 21801515

For any further questions please don't hesitate to contact us.

Sincerely,

The GPE team